

METHOD OF USAGE

- Plug the wire in to the plug set, the power indicator light flash, then into a standby situation



• diagram of plugging the power



• power switch

 Power	Power: control product power switch. When it starts, you will hear on "bi", all the function standby (power on/off, sole massage machine reset)
 Heat	Heating: heating function. Can individual control, when it works, the indicator light is blue.
 Intensity	Strength: control product massage strength, the button is used after massage function start. Three strength: green is weak, blue is medium, red is strong. The strength is default in medium when turn on product.
 Foot	Sole: sole button controls sole massage, massage focus on sole. The light is blue when it works.
 Leg	Leg: leg button controls leg massage, massage focus on leg, the light is blue when it works.
 Revive	Beautify leg: it is automatically massage leg, massage focus on leg. The light is blue when it works.
 Reflexology	Soft: it is automatically massage function. Massage focus on completely leg and sole. The light is blue when it works.



Diagram Of Turn Off The Power

Auto shut-off feature after the massager has been in use for 20 minutes. The power light is flashing and enter standby mode.

Note

1. Please switch off power and unplug power plug if not use frequently .
2. No standing during massage and before or after using to avoid damage to the machine.

PRODUCT SPECIFICATION

Name	OrthoLite
Model	RCF-O801
Rated Voltage	110-120V~ 60Hz <input type="checkbox"/> 220-240V~ 50Hz/60Hz <input type="checkbox"/>
Rated Power	100W
Rated Working Time	20Minutes
Structure Safety	<input type="checkbox"/>


SAVE THESE INSTRUCTIONS

robotouch®

OrthoLite Two-Tier

LEG BEAUTICIAN WITH ROLLER

operation instructions

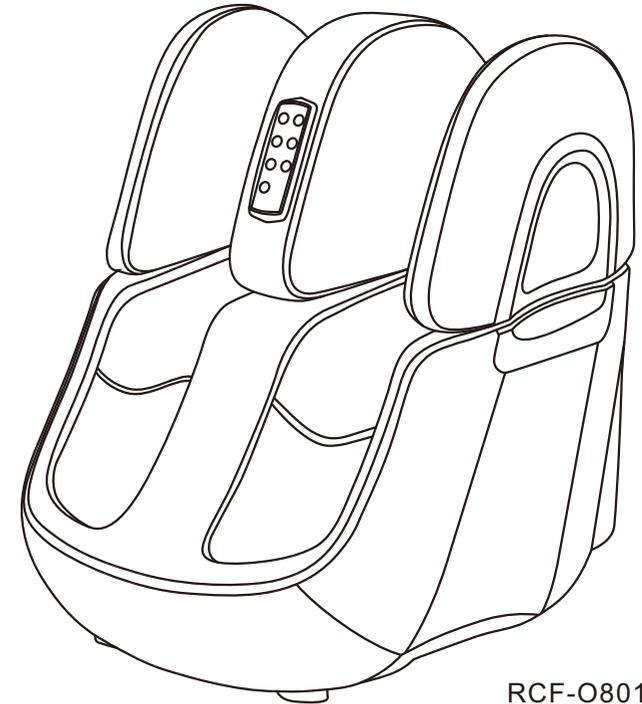

RCF-O801

Note: our company reserved the rights of product's design and explanation. Please check with us before order, we may change outlook design without notice.

IMPORTANT SAFETY WARNING



1. SafetyWarning

- Please use this product according to this Manual
- Please ensure that the mains power is suitable for this product. do not use this product if its casing or mains cable is damaged
- Do not allow any foreign matter to enter the clearance between the massagehead and the casing
- Do not allow the product to come into contact with water so as to avoid anypossible accidents that may result from a short circuit
- Please disconnect the product from the mains supply after use or before cleaningto avoid damage to the equipment or personal injury
- Do not plug or unplug this product when your hands are wet
- Do not pull the cable when plugging orunplugging, and not available strong force
- Do not damage the cable or change the circuit of the product
- Please cease using this product in case of a power cut in order to avoid anybody injury upon sudden power resumption
- Do not use any spare parts and attachments other than those recommended
- Do not use the appliance in outdoor
- Do not use this product with too much force so as to avoid injury
- Do not use this product within an hour after a meal
- Please stop using this product and consult your local supplier if you noticeanything abnormal with the product during use
- Please stop using this product and consult your doctor for adviceif you noticeanything abnormal with the product during use.

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.Children should be supervised to ensure that they do not play with the appliance.

If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.

2. Working Environment

- Do not use this product in a highly humid environment such as a bathroom
- Do not use this product immediately after there has been a sharp change in theenvironmental temperature
- Do not use this product in an environment where there are lots of dust orcorrosive gases

3. People That Should Not Use This Product

- Women that are pregnant or menstruating
- People with contagious skin conditions
- Patients with osteoporosis
- Patients with heart diseases and with electronic medical instruments like a heartpacemaker
- People in a high fever


IMPORTANT SAFETY WARNING




- Those who are injured or suffering from anybody surface disease
- Children under 14-year old or people mentally unbalanced without being Supervisedare not allowed to use this machine
- Those who need to rest according to a doctor's advice or those who do not feelwell
- Those who are wet
- Those are under the influence of alcohol or feeling unwell
- Please according to your doctor's advice if you have had an operation recently
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

FUNCTION

- 1.The new full-airbag massage leg machine has double-cut design, with a cut-over function on the top, a novel and stylish look
- 2.Airbag massage and sole roller massage: using massage chair air pressure technology, with softness strength and stereoscopic feeling .The massage techniques are flexible and changeable. The strength is from light to heavy and In-depth stimulation of acupuncture points. Sole roller massages the sole acupuncture points all of sole.
- 3.Three intensity adjustments, five kinds of programs and 17 different massage combinations: multiple program combinations for customer to choose their favor by themselves.
- 4.Multi-angle massage :the angle is adjustable (0-90°) to achieve the best results The bottom with manual flip support angle of 20 degrees, lift up to facilitate the massage of the knee.
- 5.Far infrared heating: heating function with thermal therapy effect.
- 6.Full massage: the massage range is wide, from the foot to the calf, to the knee thigh.
- 7.Safety: if the power is suddenly cut off in operation, it can be automatically deflated will not hold.
- 8.After 20 minutes, the "auto off" function will work automatically.
- 9.The kneading cloth cover with zip can be took off, and can be cleaned.

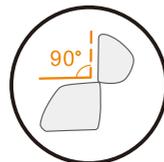

Calf turnover

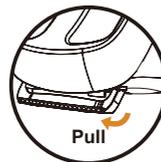


Overturn frame




Back Wheel

COMPONENT NAME AND FUNCTION

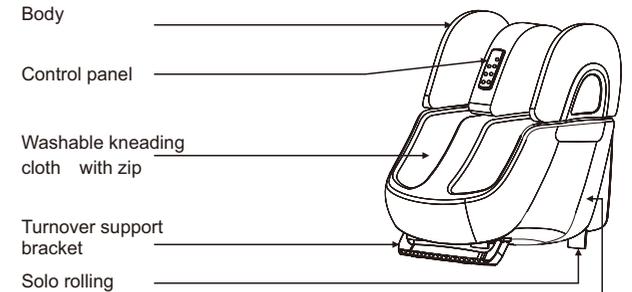

Cord and power switch

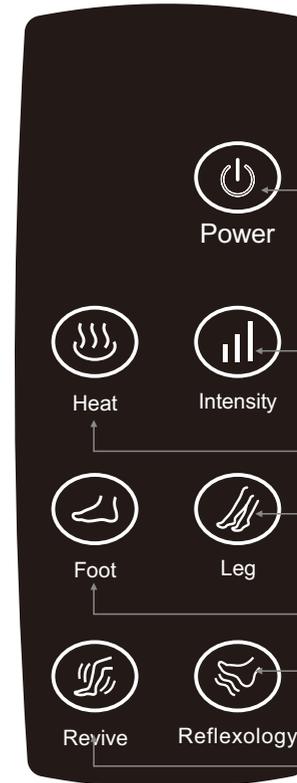

When you turn on or turn off massage function, the machine will reset.

Massage strength(Weak, medium and strong)

On/ off heating function

On/off calf massage function, massage focus on calf

On/off sole massage function, massage focus on sole.

Automatic massage function: full massage on leg and sole.

Automatic massage function: massage focus on leg